


NIR streak camera C11293-02

For picosecond time-resolved measurements of low level light emissions in the near-infrared region (1000 nm to 1650 nm) with semiconductor photocathode (InP/InGaAs)

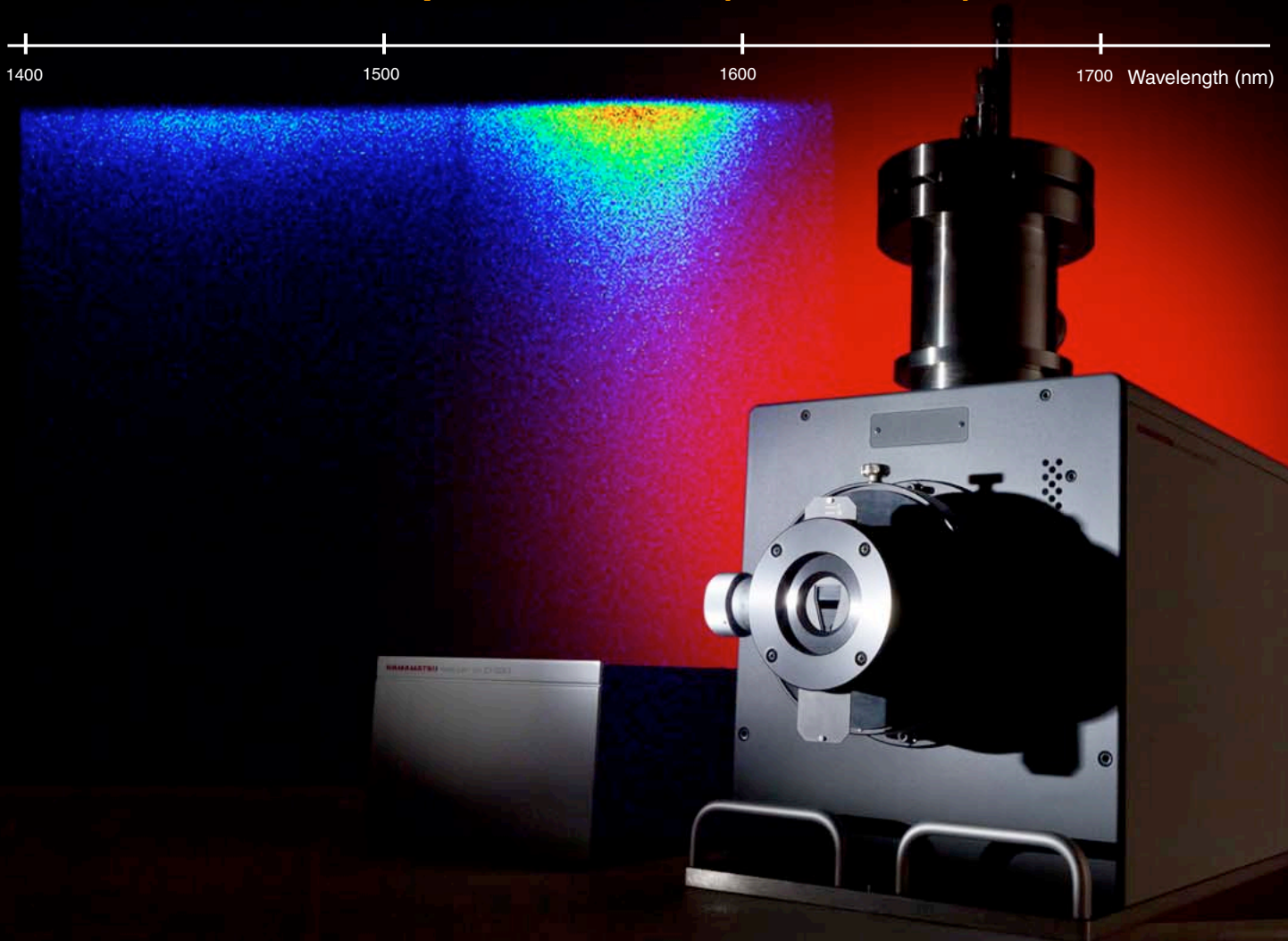

The NIR streak camera C11293-02 delivers much higher sensitivity in the near-infrared region than previously available. The C11293 utilizes a Hamamatsu streak tube with semiconductor photocathode (InP/InGaAs) to make picosecond time-resolved measurements of low-level light emissions in the near-infrared region. The C11293-02 enables simple and efficient measurements with sweep repetition rates up to 20 MHz. The streak tube provides high signal-to-noise ratio with low dark current by means of photocathode cooling to $-80\text{ }^{\circ}\text{C}$. Applications include time-resolved spectroscopy measurements of various semiconductor materials as well as time response of various light sources and devices.

Features


- Measurement wavelength: 1000 nm to 1650 nm
- Temporal resolution: 20 ps
- Multichannel measurement

Applications

- Quantum nanostructures
- Photonic crystal research
- Carbon nanotube research
- Material for photovoltaics
- Optical communication

HAMAMATSU
PHOTON IS OUR BUSINESS

System configuration example

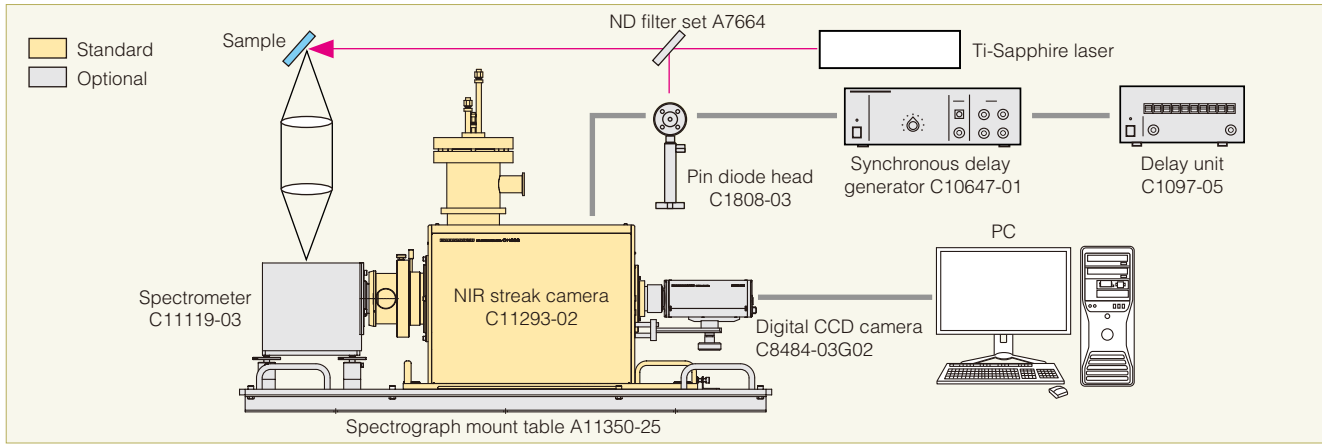


Measurement examples


Emission lifetime of quantum well (MQW488Bflat) captured with streak camera

Excitation light source: Ti-Sapphire Laser
Wavelength: 923 nm, Pulse width: 2 ps, Repetition rate: 2 MHz

Streak image

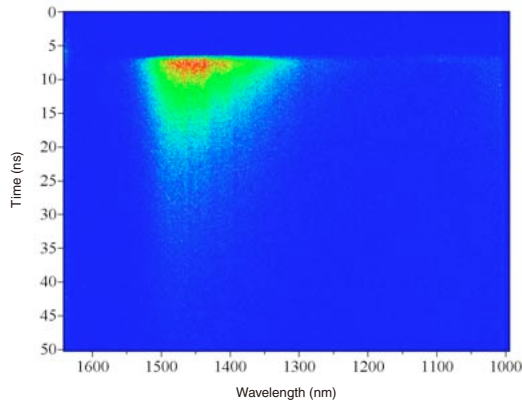

Fluorescence lifetime measurement

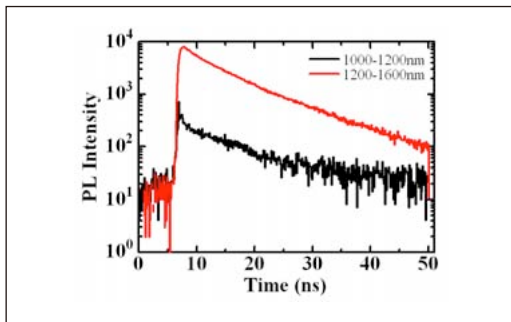


Data courtesy of Kyoto University,
Department of Electronic and Engineering, Noda Laboratory

Specifications



Type number	C11293-02
Measurement wavelength range	1000 nm to 1650 nm
Sweep time	1 ns, 2 ns, 5 ns, 10 ns, 20 ns, 50 ns, 100 ns, 200 ns 500 ns, 1 μs, 2 μs, 5 μs, 10 μs, 20 μs, 50 μs 100 μs, 200 μs, 500 μs, 1 ms, 2 ms, 5 ms, 10 ms
Temporal resolution	20 ps (FWHM)
Sweep repetition frequency (Max.)	20 MHz (1 ns, 2 ns)
Photocathode cooled temperature	-80 °C or less

Spectral response characteristics

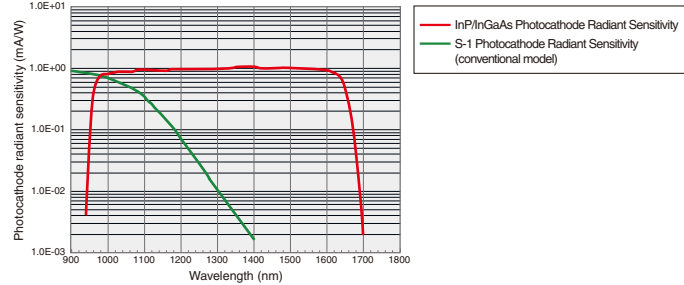

Dimensional outlines (Unit: mm)

NIR streak camera (Including spectrometer, mount table and camera)

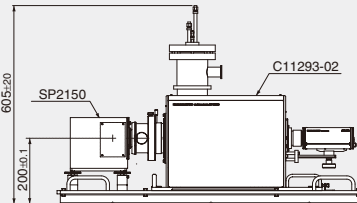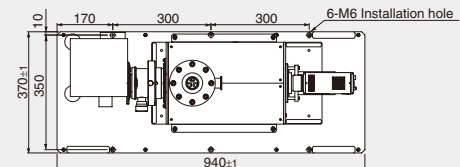

*The C11293-02 unit has a mass of approximately 29 kg.

- ★ Product and software package names noted in this documentation are trademarks or registered trademarks of their respective manufacturers.
 - Subject to local technical requirements and regulations, availability of products included in this promotional material may vary. Please consult your local sales representative.
 - Information furnished by HAMAMATSU is believed to be reliable. However, no responsibility is assumed for possible inaccuracies or omissions.
- Specifications and external appearance are subject to change without notice.
© 2015 Hamamatsu Photonics K.K.

HAMAMATSU PHOTONICS K.K. www.hamamatsu.com

HAMAMATSU PHOTONICS K.K., Systems Division

812 Joko-cho, Higashi-ku, Hamamatsu City, 431-3196, Japan, Telephone: (81)53-431-0124, Fax: (81)53-435-1574, E-mail: export@sys.hpk.co.jp

U.S.A.: Hamamatsu Corporation, 360 Foothill Road, Bridgewater, N.J 08807, U.S.A., Telephone: (1)908-231-0960, Fax: (1)908-231-1218 E-mail: usa@hamamatsu.com

Germany: Hamamatsu Photonics Deutschland GmbH, Arzbergerstr. 10, D-82211 Herrsching am Ammersee, Germany, Telephone: (49)8152-375-0, Fax: (49)8152-265-8 E-mail: info@hamamatsu.de

France: Hamamatsu Photonics France S.A.R.L.: 19, Rue du Saule Trapu, Parc du Moulin de Massy, 91882 Massy Cedex, France, Telephone: (33)1 69 53 71 00, Fax: (33)1 69 53 71 10 E-mail: infos@hamamatsu.fr

United Kingdom: Hamamatsu Photonics UK Limited, 2 Howard Court, 10 Twin Road, Welwyn Garden City, Hertfordshire AL7 1BW, UK, Telephone: (44)1707-294888, Fax: (44)1707-325777 E-mail: info@hamamatsu.co.uk

North Europe: Hamamatsu Photonics Norden AB, Torshamnsgatan 35 16440 Kista, Sweden, Telephone: (46)8-509-031-00, Fax: (46)8-509-031-01 E-mail: info@hamamatsu.se

Italy: Hamamatsu Photonics Italia S.r.l.: Strada della Moia, 1 int. 6 20020 Arese (Milano), Italy, Telephone: (39)02-93581733, Fax: (39)02-93581741 E-mail: info@hamamatsu.it

China: Hamamatsu Photonics (China) Co., Ltd.: B1201 Jiaming Center, No.27 Dongsanhuan Bellu, Chaoyang District, Beijing 100020, China, Telephone: (86)10-6586-6006, Fax: (86)10-6586-2866 E-mail: hpc@hamamatsu.com.cn

Cat. No. SHSS0014E04
JUL/2015 HPK
Created in Japan